



Endoscopic Instruments

DISPOSABLE HYBRID: D-SHAFT SERIES

The BR Surgical D-SHAFT SERIES consists of a *reusable handle* and a *disposable sheath/tip* assembly, delivering best price-quality ratio and therefore making it the most visionary concept of laparoscopic instrumentation available.

- No cleaning and sterilization of sheath is needed. All parts with imminent exposure to the patient are single use
- All dangers associated with degradation of insulating material is reduced and reliable performance ensured. Sheath/tip assembly is economically replaced after every use
- No repair or maintenance costs
- Sheath and tip are fully insulated to prevent arching
- Easy assembly/disassembly due to ball-and-socket connection between sheath and handle
- Featuring a large variety of jaws such as double action scissors, graspers and forceps

D-SHAFT SERIES



(Shafts available in 330 and 450mm lengths, in quantities of 10/box)

BR933-07410 Metzenbaum Scissor, 5mm, curved, disposable shaft



BR933-07420 Maryland Forcecp, 5mm, disposable shaft



BR933-07430 Grasping Forcecp, 5mm, atraumatic, disposable shaft



BR933-07440 Wave Grasping Forcecp, 5mm, atraumatic serrations, disposable shaft



BR932-07200
Carbon Fiber Handle
reusable



BR932-07300
Carbon Fiber Handle
with ratchet, reusable




OVER 20,000 SURGICAL INSTRUMENTS & PROCEDURE SETS • ENDOSCOPIC VIDEO EQUIPMENT & RELATED MEDICAL DEVICES

BR SURGICAL • 3500 Beachwood Court • Suite 107 • Jacksonville, FL 32224 • customerservice@brsurgical.com
Office: 904.642.1366 • Toll Free: 888.642.1366 • Fax: 904.642.1368 • www.BRSurgical.com



The BR Surgical D-SHAFT SERIES Assembly:

1. Spread open handle to extend socket
2. Slide ball of shaft into the notched area of the socket
3. Close handle and tighten rotation knob

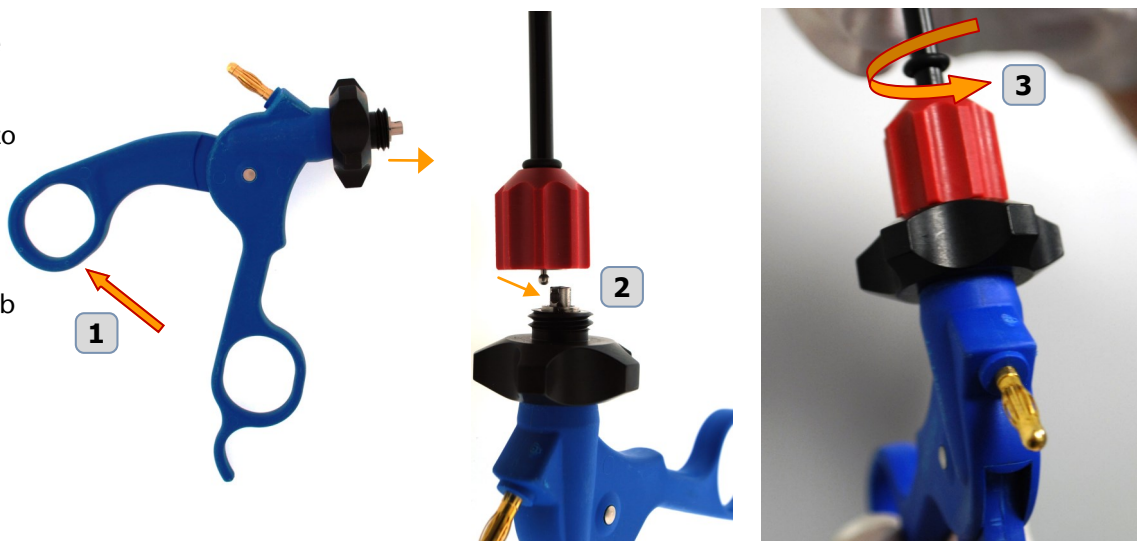

Also available: Sterilization Tray with silicon mat and instrument cradles. Call for details.

