



# Endoscopic Instruments

FULLY REUSABLE: **RS SERIES**

The BR Surgical RS SERIES is our economic line of *reusable endoscopic instruments*, consisting of a carbon fiber handle, insulated outer sheath, and a range of jaw inserts. Featuring a three-piece reusable system, the RS SERIES offers an excellent value for the money without compromises in quality.

- Fully rotatable
- Finger tang handle
- HF port
- Fully reusable
- Featuring over 130 patterns of jaws such as double action scissors, graspers and forceps

## RS SERIES (Reusable Forceps 330mm)




**BR932-07610** Metzenbaum Scissor Complete, includes handle, insulated shaft, and insert, curved, 5mm

**BR932-07620** Maryland Forcep Complete, includes handle, insulated shaft, and insert, 5mm

**BR932-07630** Grasping Forcep Complete, includes handle, insulated shaft, and insert, extremely atraumatic, serrated, 5mm

**BR932-07640** Wave Grasping Forcep Complete, includes handle, insulated shaft, and insert, with atraumatic wave-like serrations, 5mm

## REPLACEMENT INSERTS

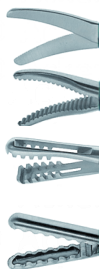

**BR932-07510-33** Metzenbaum Scissor Insert, curved

**BR932-07520-33** Maryland Forcep Insert

**BR932-07530-33** Grasping Forcep Insert, extremely atraumatic, serrated

**BR932-07540-33** Wave Grasping Forcep Insert, with atraumatic wave-like serrations



## REPLACEMENT SHEATH

**BR932-07500** Outer Sheath, insulated



**BR932-07200**  
Carbon Fiber Handle  
reusable



**BR932-07300**  
Carbon Fiber Handle  
with ratchet, reusable




OVER 20,000 SURGICAL INSTRUMENTS & PROCEDURE SETS • ENDOSCOPIC VIDEO EQUIPMENT & RELATED MEDICAL DEVICES

**BR SURGICAL** • 3500 Beachwood Court • Suite 107 • Jacksonville, FL 32224 • [customerservice@brsurgical.com](mailto:customerservice@brsurgical.com)  
Office: 904.642.1366 • Toll Free: 888.642.1366 • Fax: 904.642.1368 • [www.BRSurgical.com](http://www.BRSurgical.com)



## The BR Surgical RS SERIES Assembly:

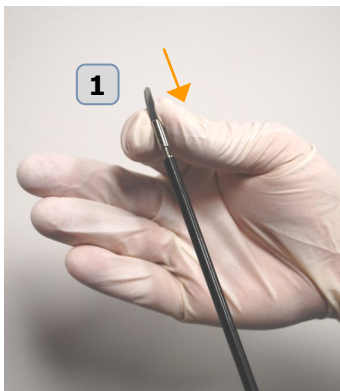

1  
Slide insert into outer sheath, thread to secure



2  
Spread open handle to extend socket



3  
Slide ball of shaft into the notched area of the socket



4  
Close handle and tighten rotation knob

Also available: Sterilization Tray with silicon mat and instrument cradles. Call for details.

